

Pets' age is different from human's!

Dog age 🐕	Human age 🧑	Cat age 🐈	Human age 🧑
1-year old	= 18-year old	1-year old	= 21-year old
3-year old	= 28-year old	3-year old	= 33-year old
7-year old	= 44-year old	7-year old	= 55-year old


- Pets' age is easy to be ignored.
- With pets' age growing the risk of acquiring diseases is increasing, especially kidney damage.
- Regular urine test is essential to the aging pets.

The significance of urine test 🐾

Urine is generated by the kidney, which contains different kinds of metabolites. Some diseases could cause the level change of the metabolites, especially kidney damage. Urine test is to examine the levels of these metabolites in order to help veterinarians to diagnostic the diseases.

The simplest way is to observe the colors of the urine. The normal color is light yellow, the obvious color change need to draw the attention to pets owners, which probably reflect the relative diseases. The best way is to bring the pets to the clinic for advices.

Generally veterinarians would recommend to make the urine test for better diagnostics!

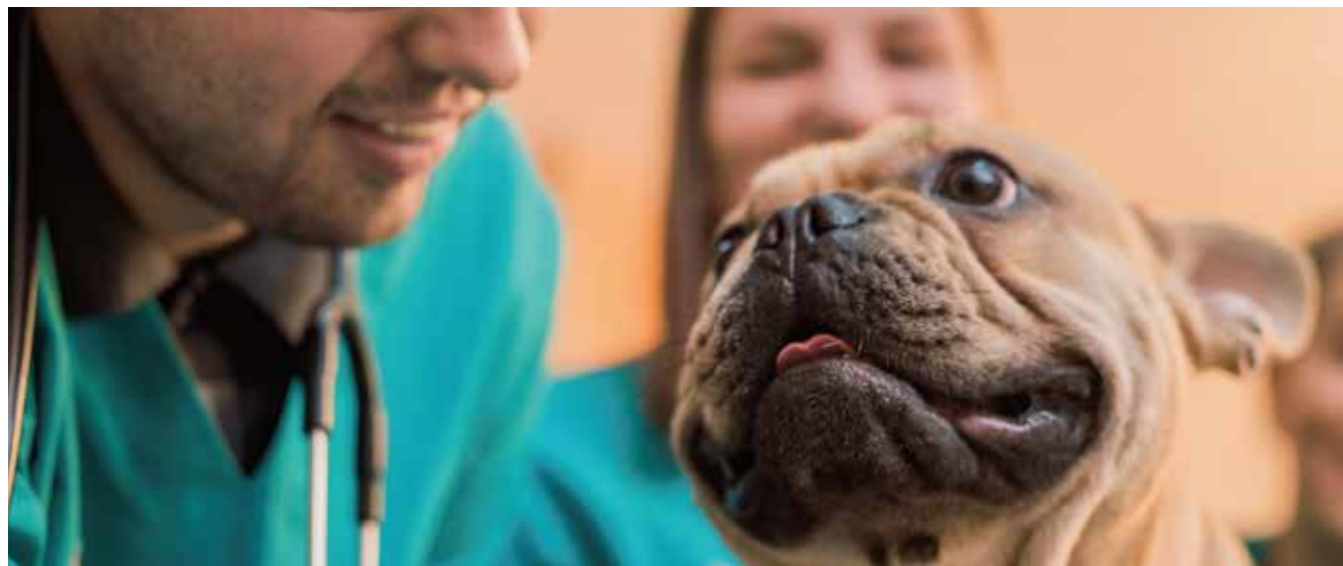

UA-60V Veterinary Urine Analyzer 🐾



Portable design
compact size



More parameters
two strips



Easy operation
efficient workflow

		Parameters
Strips-15	Strips-11	LEU
		URO
		BIL
		GLU
		ASC
		SG
		KET
		NIT
		pH
		BLO
		PRO
	CRE	
	MA	
	CA	
	ACR	



Basic feature

- Throughput: ≥120 T/H (Fast mode); ≥40 T/H (Normal mode)
- Storage: 2,000 results
- CV≤1.0%

Technical specifications

- Principle: Reflectance measurement of cold light source
- Wavelength: 525nm, 635nm

Parameter

- LEU , URO , BIL , GLU , ASC , SG , KET , NIT , pH , BLO , PRO , CRE , MA , CA , ACR

Species

- Cat, Dog

Printer

- Thermal printer

Operating environment

- Temperature: 15-30℃
- Humidity: 30%-85%

Weight and dimensions

- Weight: 0.2kg
- Dimensions: 120mm X 80mm X 29.5mm

Strips

- Type: Strips-11 and Strips-15
- Volume: 50 pcs/pack
- Shelf life: 24 months
- Open bottle stability: 3 months
- Storage conditions: Store in a dry place at 0-30℃ and avoid exposure to direct sunlight

UA-60V

Veterinary Urine Analyzer

Kidney Care Within Reach

